Lace Align[™] Spinal Orthoses

Keeping you in your comfort zone.

Spinal Orthoses

Product Number	S/M	L/XL
Lumbar Orthosis (LO)	35-3111	35-3113
Suggested Billing Code: L-0627 (Custom Fit), L-0642 (OTS)		

Lumbosacral Orthosis (LSO) 35-3201 35-3203 Suggested Billing Code: L-0631 (Custom Fit), L-0648 (OTS)

Lumbosacral Plus Orthosis (LSO PLus) 35-3211 35-3213 Suggested Billing Code: L-0637 (Custom Fit), L-0650 (OTS)

 Sizing Guide
 S/M
 L/XL

 30" - 45"
 40" - 60"

Measurement taken from waist or hip circumference. Use the larger of the two circumferences when selecting brace size.

Accessories

A. Replacement Back Panel (LO Model)	35-3901
B. Replacement Back Panel (LSO or LSO Plus Model)	35-3902
C. Replacement Side Panels, pair (LSO Plus Model)	35-3903
D. Lace Align Extension Kit (15")	35-3904
E. Gel Pouch with 9"x11" Gel Pack	35-3911
Gel Pouch Only	35-3910
9" x 11" Gel Pack (each)	72-7006
9" x 11" Gel Pack (box of 10)	72-7006-010





CORFLEX.

Designed for life.

900 426 7252

CORFLEX.COM

669 EAST INDUSTRIAL PARK DRIVE / MANCHESTER, NH 03109 / FAX 603 623 4111 / EMAIL SERVICE@CORFLEX.COM



Lace Align™ Spinal Orthoses

Keeping you in your comfort zone.






Introducing Outlast® proactive heat management technology

CORFLEX®

Designed for life.

Power Pull™ System

The key to an effective back brace is delivering support to the spinal column. Our Power Pull System allows patients to achieve this with minimal effort, stabilizing internal structures and helping reduce intradiscal pressure through circumferential compression.

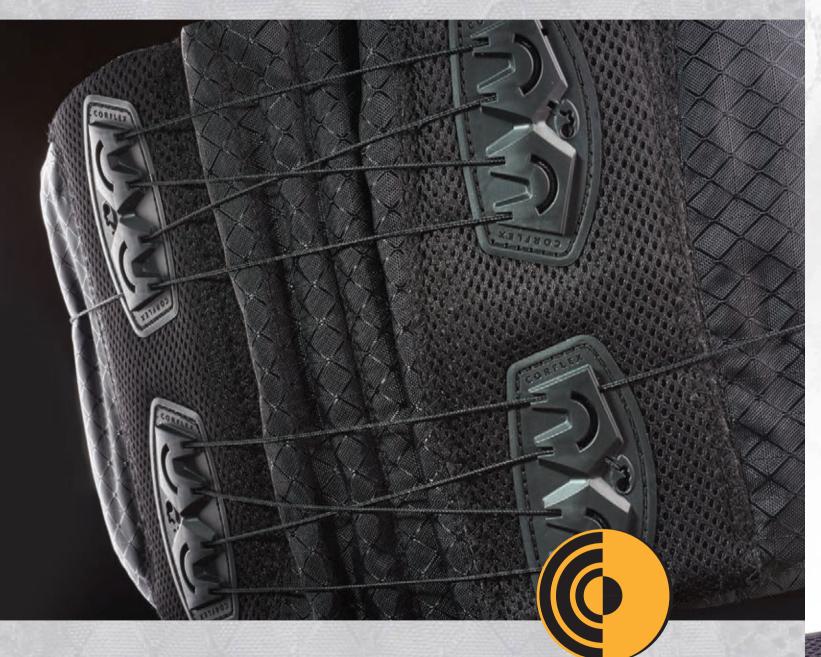

CORFLEX®

Designed for life.

Corflex® and Outlast® Keeping you in your comfort zone.

How does Outlast Work?

Outlast technology delays the onset of sweat, proactively buffering temperature swings to help prevent perspiration. This patented phase-change material absorbs, stores and releases excess body heat to maintain balance, helping the body minimize temperature swings to conserve metabolic energy.


Creating the ideal comfort zone...

We've coupled the proactive benefits of Outlast with our reactive CoolMeshTM fabric to deliver the ultimate in orthopedic brace comfort. Outlast starts working the minute you apply the brace, keeping you cool and dry. As your activity level increases, our breathable CoolMesh steps in to help wick moisture and encourage airflow, creating the ideal comfort zone.

An ideal skin temperature is maintained when the microclimate between your skin and the brace is balanced

Excess heat generated by the body is absorbed into the Outlast® Thermocules™

Breathable CoolMesh™ helps wick moisture and encourage airflow

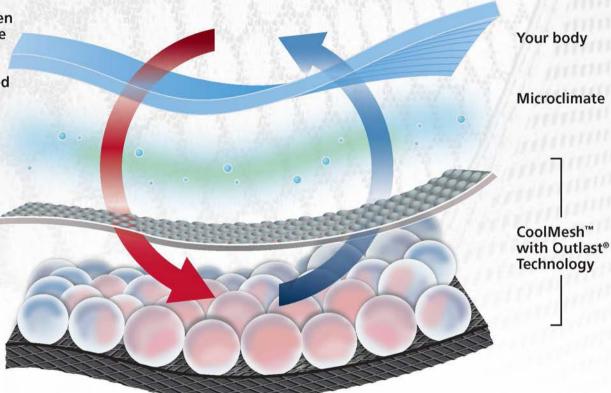



Lace Align[™] Spinal Orthoses

Rigid Back Panel

step-down modularity


a heat gun for added customization

and can be shaped with

Removes easily for

> Designed to stabilize and decompress the spine in patients recovering from > Indications: lower back pain, sprains or strains, spinal stenosis, surgery or injury, this innovative brace incorporates Outlast[®], a hi-tech fabric that provides proactive temperature regulation, managing heat and moisture for all-day comfort. Removable posterior panel allows "step-down" to a traditional back brace, offering moderate support and increased patient value. Also available is an optional Gel Pouch that attaches easily to deliver cold or hot therapy following surgery, injection or injury.

spondylosis, spondylolysis, spondylolisthesis, degenerative disc disease, bulging or herniated discs, burst fractures, facet syndrome and postoperative stabilization following laminectomy, foraminotomy, IDET procedures, discectomy, fusion or decompression.

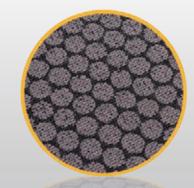

CoolMesh with

Outlast

Proactively balances temperature swings to help prevent sweating and improve patient compliance

CoolMesh Fabrile

Breathable mesh helps wick moisture and encourage airflow keeping you cool and dry

Sizing Guides

Embroidered sizing marks take the guesswork out of belt length adjustment



Lumbosacral Plus Orthosis (LSO Plus) L0637/L0650

Delivers motion control in multiple planes with the addition of movable lateral side panels



Universal Sizing

Two sizes fit waist circumference ranges of 30"-60" using the Lift-and-Slide adjustment – no trimming required



Power Pull Glosures

Durable molded pull tabs allow the patient to effortlessly apply and adjust the desired level of compression



Hand Loops

Located on both front panels to promote independent brace application and removal when grip strength is impaired



Included in the LSO Plus model for added lateral stabilization



Shaped and slotted to help eliminate irritation when sitting or bending



CORFLEX_®

Designed for life.

Patent pending.