

Motion OA 2

Osteoarthritis Knee Brace



CORFLEX®

Designed for life.

CORFLEXGLOBAL.COM


| Easy to Apply, Stock and Wear.

The Motion OA 2 brace is a lightweight single-hinge brace that is ideal for patients with mild to moderate symptoms. It provides increased rehabilitative and functional support for diagnoses that require reduced medial or lateral joint loads. The Motion OA 2 is laterally applied for both a medial and lateral OA, making it ideal for bilateral brace users.

The adjustable Q-hinge design builds correction into the brace to relieve the pressure and pain of knee osteoarthritis. Easy to apply and breathable, the high stretch fabric helps maximize fit, compression and suspension.

ADDITIONAL FEATURES:

- Universal sizing promotes inventory reduction - only three sizes to fit all of your patients; quick release buckles allow for easy application and removal; laterally applied brace helps unload the knee without interfering with gait
- Indicated for mild to moderate unicompartmental osteoarthritis (OA); Osteochondritis Dissecans (OCD); meniscus tears
- Constructed of non-neoprene, breathable fabric for improved comfort; lycra popliteal; padded hinge and 3/8" thick condyle pad

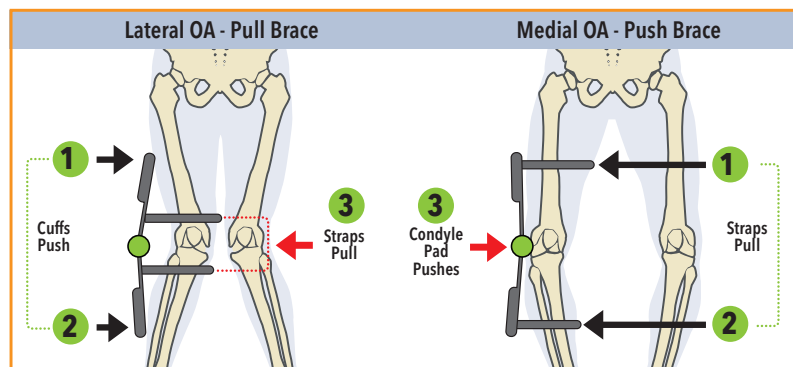

Q-Hinge angle adjustment delivers three-point leverage for maximum relief

The Motion OA 2 Single upright, off-loading brace helps restore comfort in your everyday life. The wrap style brace is easy to apply and provides the right fit, easily conforming to an individual's anatomy.

Motion OA 2 Osteoarthritis Knee Brace

Finger Loops

Provide ease of application.

Improved Comfort

Constructed of non-neoprene, breathable fabric with Lycra material in back of knee to reduce bunching and popliteal irritation.

Adjustable Knee Extension

Metal knee hinges feature movable stops that control knee extension from 0° to 25° in 5° increments and flexion from 45° to 90° in 5° increments.

Push Buttons

Uses quick release buckles for hassle-free application.

Easy to Apply


Wrap style brace features low profile, flexible cuffs that resist rotation.

Q-Hinge

Varus/valgus adjustment provides compartmental off-loading to help maximize pain relief.

Anatomically Contoured

High stretch fabric provides greater fit and compression.

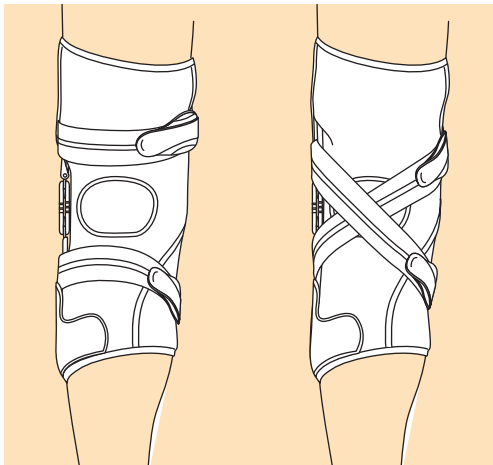

Straight Strapping
for Medial OA
Application

Criss-Cross Dynamic
Strapping for Lateral OA
Application

MOTION OA 2

Product Number

	S/M	M/L	L/XL
Right	81-0222	81-0224	81-0226
Left	81-0232	81-0234	81-0326

Sizing Guide

	S/M	M/L	L/XL
6" Above Patella	14"-18½"	18½"-25½"	25½"-32"
6" Below Patella	11"-14"	14"-19"	19"-25"

Suggested Billing Codes L-1843 (Custom Fit) or L-1851 (OTS)
Knee Orthosis Single Upright



CORFLEX GLOBAL | 669 EAST INDUSTRIAL PARK DRIVE | MANCHESTER, NH 03109
CORFLEXGLOBAL.COM | TOLL FREE: 800 426 7353 | FAX: 603 623 4111