

CoolTex™ Ag 6" Thumb Spica

Low profile design and flexible support provide secure and comfortable immobilization



CORFLEX®


Designed for life.

CORFLEXGLOBAL.COM



Our 6" CoolTex Ag Thumb Spica is designed to provide secure and comfortable immobilization of the wrist and thumb while preserving finger mobility. Single piece design minimizes seams to ensure maximum patient comfort while providing uniform compression. Wide elastic wrist strap provides additional compression that is easily controlled with a hook and loop closure. Splint features removable, moldable, plastic palmar, thumb and dorsal stays for optimal stabilization. Elastic wrist strap provides additional support and compression.

Lightweight and breathable CoolTex Ag material allows for greater comfort during extended wear, helping wick moisture away from the skin. Material also contains silver yarns that provide antimicrobial protection to the fabric to help reduce odors and protect fibers from premature breakdown.

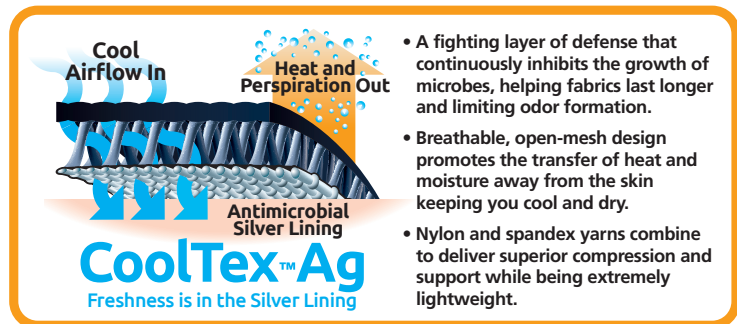

ADDITIONAL FEATURES:

- Shorter length is ideal for use during daily activities; removable stays allow for customized support of the thumb and wrist independently; elastic wrist strap allows variable compression with an easy-to-adjust hook and loop closure
- Indicated for deQuervain's Syndrome, tendinitis, thumb strains or sprains and post-cast removal
- Constructed of breathable, moisture-wicking nylon spacer material; removable, semi-rigid HPDE and LDPE dorsal stays and thumb spica and ABS palmar stay

CoolTex™ Ag 6" Thumb Spica



Dual dorsal stays made from semi-rigid LDPE and HDPE plastic are removable to allow varied support based on activity.

Product Number		Size	Measurement (Wrist Circumference)
Right	Left		
74-5640	74-5650	XS	5" - 5 3/4"
74-5641	74-5651	S	5 3/4" - 6 1/2"
74-5642	74-5652	M	6 1/2" - 7 1/4"
74-5643	74-5653	L	7 1/4" - 8"
74-5644	74-5654	XL	8" - 8 3/4"
74-5645	74-5655	2XL	8 3/4" - 9 1/2"

Suggested Billing Codes

L-3807 (Custom Fit) WHFO, no joint, prefabricated or
L-3809 (OTS) WHFO without joints prefabricated



Rigid yet formable ABS plastic palmar stay can be contoured for an individualized fit.

ALSO AVAILABLE

CoolTex™ Ag 6" Wrist Splint

Our 6" CoolTex Ag Wrist Splint is ideal for all-day wear. Support levels can be easily adjusted by removing stays or changing compression support at the wrist.

ASK ABOUT OUR FULL LINE OF
COOLTEX AG PRODUCTS





Semi-rigid LDPE plastic thumb stay can be removed and formed for customized thumb support.



CORFLEX GLOBAL | 669 EAST INDUSTRIAL PARK DRIVE | MANCHESTER, NH 03109
CORFLEXGLOBAL.COM | TOLL FREE: 800 426 7353 | FAX: 603 623 4111